


OsteoCentric®

Small Fragment Set

Technique Guide

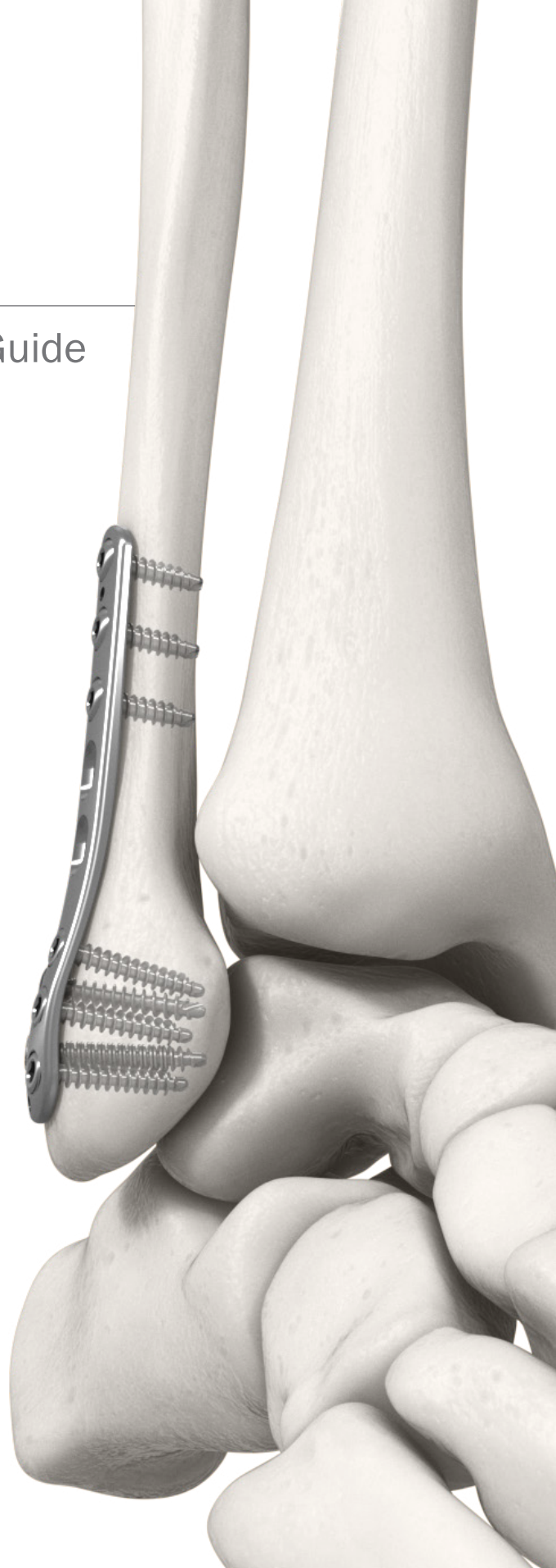


OsteoCentric®
TECHNOLOGIES

OsteoCentric® Small Fragment Set Featuring UnifiMI® Technology

Table of Contents

UnifiMI® Technology	3
Indications for Use	4
Contraindications	5
Compression, 1/3 Tubular, & Reconstruction Plate Technique	5
Lateral Distal Fibula Plate Technique	7
Lateral Locking Distal Fibula Plate Technique	7
Fastener Only Technique	8
Long Fastener Technique	8
Implants & Instruments List	9
4.0 Super Fasteners & Instruments	11

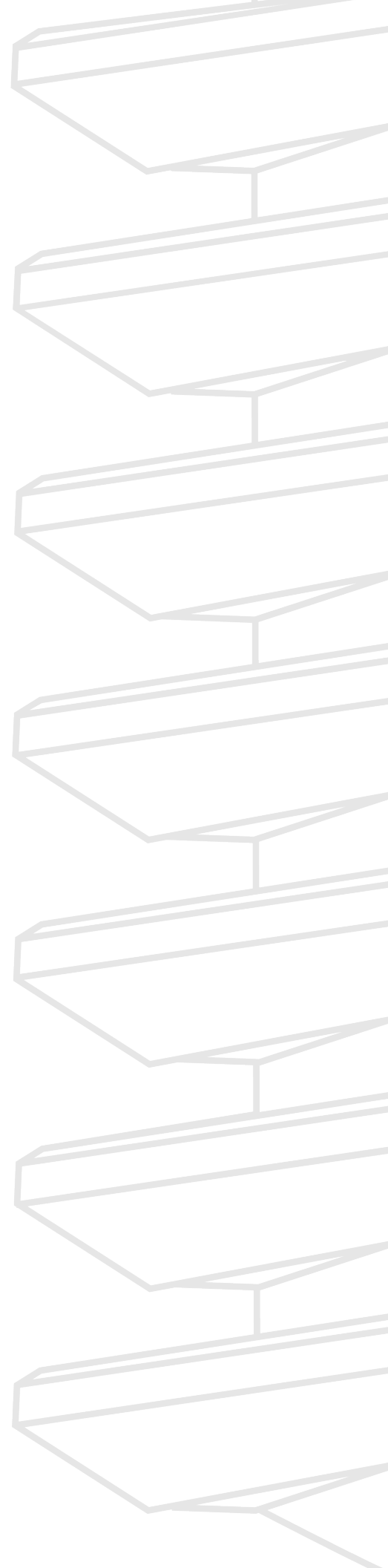



Mechanical Integration (MI) is a revolutionary new method of instantly securing implants to bone utilizing OsteoCentric's proprietary thread geometry called UnifiMI - providing superior primary stability of any implant in normal or compromised bone. Unique thread geometry instantly and circumferentially interlocks with bone by entrapping and containing bone between the thread form. This mechanical interlocking creates a structural and functional connection between an implant and bone which performs similar to Biological Integration (Osseointegration). Unlike Biological Integration (Osseointegration), the UnifiMI connection to bone is instant vs. delayed, removable vs. permanent, mechanical vs. biological.

The screw fasteners in the Small Fragment Set are designed to preserve the integrity of and instantly interlock to the bone with its unique thread geometries. (Figure 1)

Thin Profile, Interlocking Strength

UnifiMI's interlocking threads are designed to increase construct stability while maintaining the benefits of low-profile plates.

Versatile Fixation

The set includes Compression, Reconstruction, 1/3 Tubular, and Distal Fibula plates along with 3.5mm, 3.3mm, and 2.7mm interlocking screw fasteners that work in any plate hole for customized constructs for each patient.


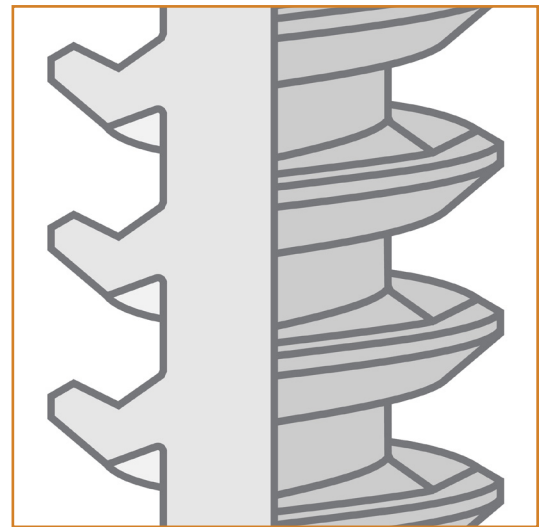


Figure 1

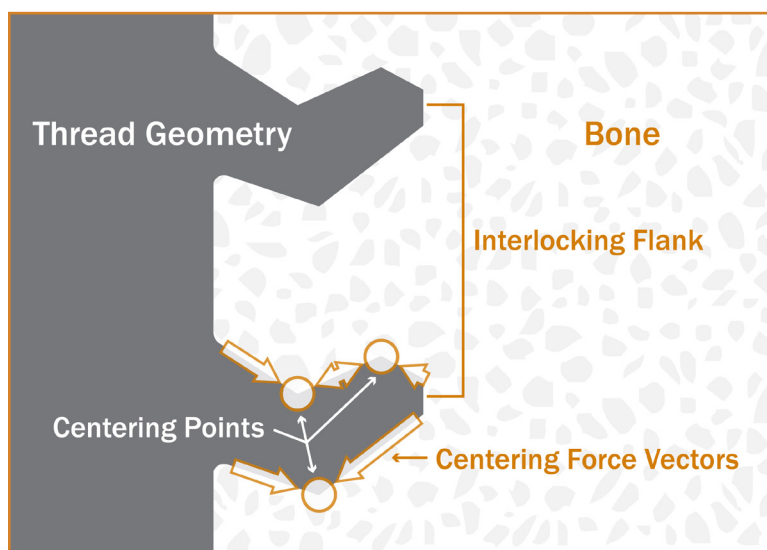

UnifiMI Technology

The unique UnifiMI technology provides the benefits of robust fixation and immediate mechanical interlocking with the bone without the need for special techniques or instruments.

The 3.5mm Screw Fastener provides 80° of total angulation proximal-distal and 14° of total angulation axially to allow for variation in the fixation needs. (Figures 2 and 3)


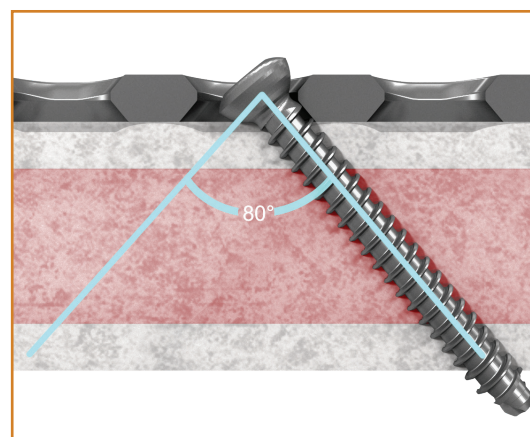


Figure 2

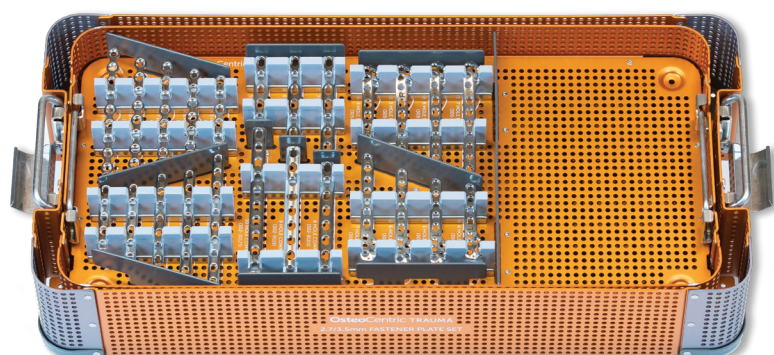

Fastener Plate Set

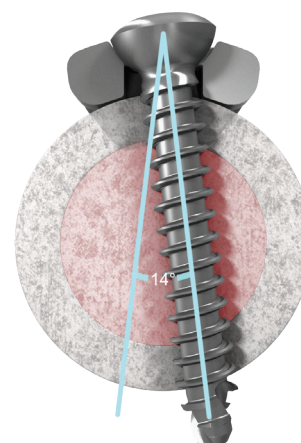


Figure 3

The OsteoCentric Small Fragment set uses all standard instrumentation familiar to orthopedic surgeons. Fasteners and implants are currently available in stainless steel.

Indications For Use

The OsteoCentric Trauma 2.7mm Screw Fastener is intended for fractures and osteotomies of small bones and bone fragments, including the foot, ankle, and hand in adults and in both children (2-12 years) and adolescents (12-21 years) in which growth plates have fused or in which growth plates will not be crossed by screw fixation.

The OsteoCentric Trauma 3.5mm Screw Fasteners are intended for fixation of fractures, osteotomies and non-unions of the clavicle, scapula, olecranon, humerus, radius, ulna, pelvis, tibia, calcaneus, femur and fibula in adults and in both children (2-12 years) and adolescents (12- 21 years) in which growth plates have fused or in which growth plates will not be crossed by fastener fixation.

The OsteoCentric Trauma Bone Plate and Screw System is intended for fixation of fractures, osteotomies, and non-unions of the clavicle, scapula, olecranon, humerus, radius, ulna, pelvis, tibia, and fibula, including, periarticular and intraarticular fractures.

Contraindications

The physician's education, training, and professional judgment are necessary to determine the appropriate treatment protocol and patient selection. Contraindications may be relative to each patient, and clinicians should always consider all risks and possible reactions when considering the proper treatment protocol.

Specific contraindications include:

- Active or latent infection
- Insufficient quantity or quality of bone/soft tissue
- Material sensitivity – If suspected, tests should be performed prior to implantation.
- Sepsis
- Patients who are unwilling or incapable of following postoperative care instructions.
- Spinal fixation – This device is not intended for screw attachment or fixation to the posterior elements (pedicles) of the cervical, thoracic, or lumbar spine.


Compression, 1/3 Tubular, & Reconstruction Plate Technique

Select the type of plate appropriate for the fracture and anatomic area.

If bending is needed, use the Bending Iron L and R or the Bending Iron Reconstruction to customize the plate contour.

Hold the plate to the bone using the Reduction Forceps provided.

Instruments

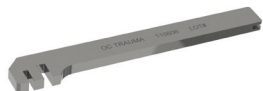

Bending Iron
Left: 110605
Right: 110606



Bending Iron Reconstruction
110610



Reduction Forceps, Points
Narrow: 110630
Broad: 110631



Reduction Forceps, Serrated
110632

Drill for the first screw fastener in the neutral position (*Figure 4*) with the appropriate drill bit for the screw size selected.

Screw Fastener

2.7mm
3.5mm

Drill Bit

2.0mm
2.5mm

Measure using the Depth Gauge. (*Figure 5*)

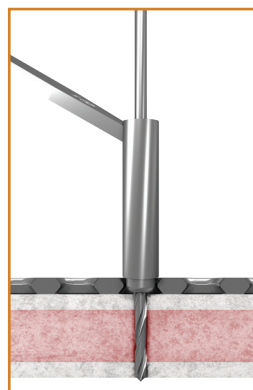


Figure 4


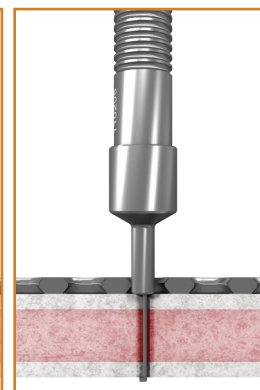



Figure 5

Instruments

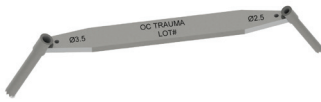

Drill Sleeve

3.5mm/2.5mm: 110520
2.7mm/2.0mm: 110521

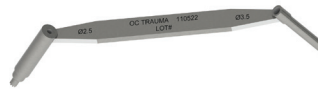

Universal Drill Sleeve

3.5mm/2.5mm: 110522



Insert Drill Sleeve

3.5mm/2.5mm: 110524



Drill Bit

2.0mm, 100mm: 110007
2.5mm, 110mm: 110013
2.5mm, 145mm: 110014



Depth Gauge

60mm: 110205

1. If dense, hard bone is detected, please check to be sure the proper drill bit size has been utilized. If it has, it is recommended to use the appropriately sized tap to prepare the bone for screw insertion.
2. Insert the screw fastener using the preferred 2.5mm hex screwdriver.
3. Drill for a screw fastener on the opposite side of the fracture. For compression at the fracture with Compression Plates, drill eccentrically (Figure 6) with the drill bit away from the fracture.
4. Measure using the Depth Gauge and insert the screw fastener using the preferred 2.5mm hex screwdriver.
5. For additional compression at the fracture site, lag screw technique can be used through the plate by overdrilling the near cortex (Figure 7) with the 3.5mm Drill Bit or 2.7mm Drill Bit.
6. Insert remaining screw fasteners as appropriate for the fracture.



Figure 6



Figure 7

Instruments

Drill Bit

2.7mm, 100mm: 110016
3.5mm, 110mm: 110028



Lateral Distal Fibula Plate Technique

The Distal Fibula Plate is implanted using the previous techniques. (Figure 8) Every screw hole can accept any screw size, however, 2.7mm or 3.3mm Screw Fasteners are often used in the distal cluster with 3.5mm used in the shaft (Figure 9).

Screw Fastener	Drill Bit
2.7mm	2.0mm
3.3mm	2.0mm
3.5mm	2.5mm

Lateral Locking Distal Fibula Plate Technique

Leveraging the UnifitTM proprietary design, the Locking Distal Fibula Plate allows surgeons to place both locking or non-locking fasteners into the same hole, resulting in lower head prominence (Figure 10). Determine the required length of the fastener by using the scale on the 2.0mm calibrated drill guide (Figure 11).

Screw Fastener	Drill Bit
3.3mm	2.0mm



Figure 8

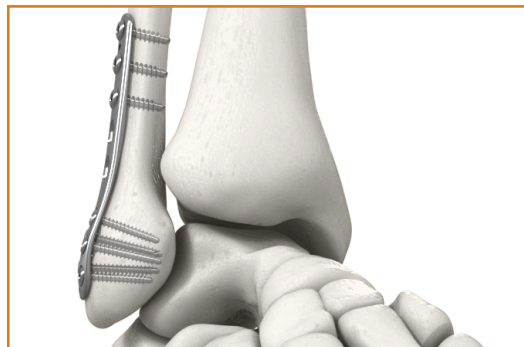


Figure 9

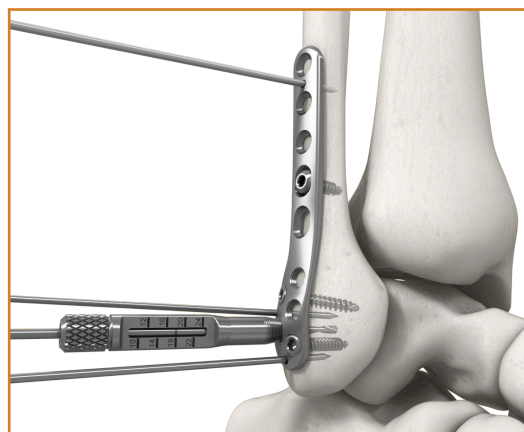


Figure 10



Figure 11

Fastener Only Technique

Screw Fastener implantation uses standard techniques with the appropriate drill bit for the implant chosen.

Screw Fastener	Underdrill	Overdrill
2.7mm	2.0mm	2.7mm
3.5mm	2.5mm	3.5mm
4.0mm	2.5mm	4.0mm

Note: 4.0mm super fasteners are optimized to perform in metaphyseal region, or in the case of weak or compromised bone.

After drilling and before measuring, the Countersink can be used to reduce the profile of the head on the near cortex.

Instruments



- Countersink: 110500

Note: If during any fastener insertion technique excessive insertion torque is encountered, it is highly recommended to remove fastener and utilize the correctly sized tap to prepare the bone for fastener insertion.



Long Fastener Technique

Fixation needs may require the use of extra-long fasteners and corresponding instruments. The following are part of the Long Fastener Instrument Set:


Instruments

- Drill Bit
 - 3.5mm Drill Bit 110mm: 110028
 - 2.5mm Drill Bit Calibrated 200mm: 110040
- Hex Driver
 - Long Hex Driver: 110104
 - Long Hex Driver Fixed 300mm: 110703
 - 2.5mm Hex Driver Sleeve: 110701
- Depth Gauge For Small Frag 150mm: 110206
- 3.5mm Cortical Screw Fastener Tap 200mm: 110454
- Screw Fastener Holding Forceps: 110705
- 2.0mm Drill 145mm, Calibrated: 110039
- 2.7mm Cortical Screw Tap 145mm: 110402
- 2.7mm Drill Bit 100mm: 110016

Note: If during any long fastener insertion technique excessive insertion torque is encountered, it is highly recommended to remove fastener and utilize the correctly sized tap to prepare the bone for fastener insertion.

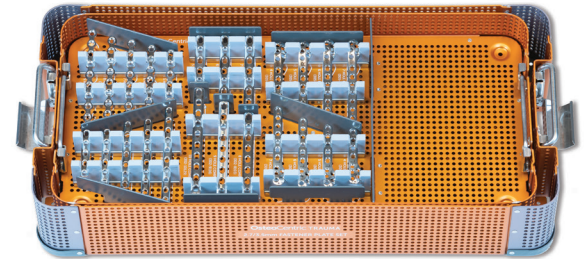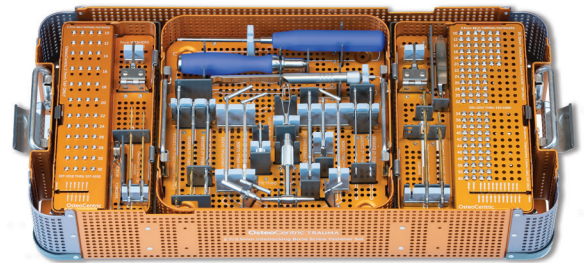

OsteoCentric Trauma Small Fragment Set

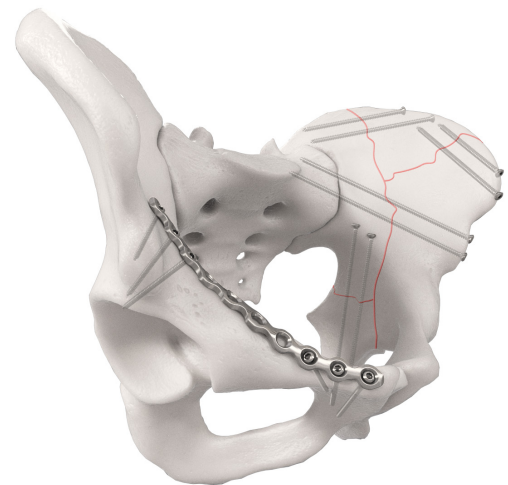

Small Fragment Set Implants

3.5mm Compression Plate

Part #	Holes	Length (mm)
302-564	5	64
302-677	6	77
302-790	7	90
302-8103	8	103
302-9116	9	116
302-10129	10	129



1/3 Tubular Plate With Collar

Part #	Holes	Length (mm)
303-557	5	57
303-669	6	69
303-781	7	81
303-893	8	93
303-10117	10	117



Straight Reconstruction Plate

Part #	Holes	Length (mm)
304-565	5	65
304-678	6	78
304-791	7	91
304-8104	8	104
304-9117	9	117



K-Wires

Part #	Diameter	Length (mm)
130005	1.25mm	150
130010	1.6mm	150
130015	2.0mm	150

2.7mm/3.5mm Lateral Distal Fibula Plate

Part# L	Part# R	Holes	Length (mm)
305-470L	305-470R	4	70
305-580L	305-580R	5	80
305-690L	305-690R	6	90
305-7100L	305-7100R	7	100



2.7mm/3.5mm Locking Lateral Distal Fibula Plate

Part# L	Part# R	Holes	Length (mm)
314-470L	314-470R	4	70
314-580L	314-580R	5	80
314-690L	314-690R	6	90
314-7100L	314-7100R	7	100
314-9120L	314-9120R	9	120
314-11140L	314-11140R	11	140



2.7mm Screw Fastener Self-Tapping

10mm - 32mm, 2mm increments

327-10xx

2.7mm Long Screw Fastener Self-Tapping

34mm - 50mm, 2mm increments

50mm - 110mm, 5mm increments

327-1xxx



3.3mm Screw Fastener Self-Tapping

12mm - 22mm, 2mm increments

333-10xx

3.3mm Locking Screw Fastener Self-Tapping

12mm - 22mm, 2mm increments

333-60xx



3.5mm Screw Fastener Self-Tapping

10mm - 50mm, 2mm increments

335-10xx

3.5mm Long Screw Fastener Self-Tapping


52mm - 88mm, 2mm increments

90mm - 130mm, 5 mm increments

335-1090-1130

7mm Washer

300-0701

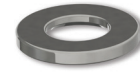

9mm Washer

300-0901

Full Instrument List

Part #	Description	Part #	Description
110007	2.0mm Drill Bit 100mm	110625	Periosteal Elevator 6mm Width
110013	2.5mm Drill Bit 110mm	110630	Reduction Forceps Narrow
110014	2.5mm Drill Bit 145mm	110631	Reduction Forceps Broad
110016	2.7mm Drill Bit 100mm	110632	Reduction Forceps Serrated
110028	3.5mm Drill Bit 110mm	110701	2.5mm Hex Driver Sleeve
110040	2.5mm Drill Bit Calibrated 200mm	110702	2.5mm Hex Driver with Silicone Handle
110045	2.0mm Drill Bit Calibrated 110mm	110703	Long Hex Driver Fixed 300mm
110103	2.5mm Hex Driver QC Shaft 110mm	110705	Screw Fastener Holding Forceps
110104	Long Hex Driver	130005	1.25mm X 150mm K-Wire Smooth
110205	Depth Gauge for Small Frag 60mm	130010	1.6mm X 150mm K-Wire Smooth
110206	Depth Gauge For Small Frag 150mm	130015	2.0mm X 150mm K-Wire Smooth
110302	AO QC Handle	A2200	Small Fragment Fastener Set Outer Case
110310	QC T-Handle	A2200-01	Lid
110422	2.7mm Screw Fastener Tap 110mm	A2200-02	Small Fragment Set Instrument Tray Upper
110424	3.5mm Screw Fastener Tap 110mm	A2200-03	Small Fragment Set Instrument Tray Lower
110444	3.5mm Screw Fastener Tap 135mm	A2200-04	3.5mm Small Fragment Fastener Caddy
110454	3.5mm Cortical Screw Fastener Tap 200mm	A2200-06	2.7/3.3mm Small Fragment Fastener Caddy
110500	3.5mm Countersink	A2200-07	2.7/3.3/3.3mm Locking Fastener Caddy
110520	3.5mm/2.5mm Drill Sleeve	A2300	Small Fragment Plate Set Outer Case
110521	2.7mm/2.0mm Drill Sleeve	A2300-01	Small Fragment Plate Set Instrument Tray
110522	3.5mm/2.5mm Universal Drill Sleeve	A2300-02	Small Fragment Plate Set Implant Tray
110524	3.5mm/2.5mm Insert Drill Sleeve	A2300-03	Locking LDFP and Instrument Tray
110535	2.0mm Threaded Drill Guide Calibrated	A2300-04	Locking LDFP 9 Hole, 11 Hole Plate Tray
110600	Sharp Hook	A2400	Long Fastener Outer Tray
110605	Bending Iron	A2400-02	Long 3.5mm Fastener Caddy
110606	Bending Iron	A2400-03	Long Fastener Instrument Tray
110610	Bending Iron Reconstruction	A2400-01	Long 2.7mm Fastener Caddy
110620	Hohmann Retractor 8mm Width	110039	2.0mm Drill 145mm, Calibrated
110621	Hohmann Retractor 15mm Width	110402	2.7mm Cortical Screw Tap 145mm
		110016	2.7mm Drill Bit 100mm

4.0mm Super Fasteners

Part Number	Description	Length
340-2012	4.0 Bone Screw Super Fastener ST	12mm
340-2014	4.0 Bone Screw Super Fastener ST	14mm
340-2016	4.0 Bone Screw Super Fastener ST	16mm
340-2018	4.0 Bone Screw Super Fastener ST	18mm
340-2020	4.0 Bone Screw Super Fastener ST	20mm
340-2022	4.0 Bone Screw Super Fastener ST	22mm
340-2024	4.0 Bone Screw Super Fastener ST	24mm
340-2026	4.0 Bone Screw Super Fastener ST	26mm
340-2028	4.0 Bone Screw Super Fastener ST	28mm
340-2030	4.0 Bone Screw Super Fastener ST	30mm
340-2032	4.0 Bone Screw Super Fastener ST	32mm
340-2034	4.0 Bone Screw Super Fastener ST	34mm
340-2036	4.0 Bone Screw Super Fastener ST	36mm
340-2038	4.0 Bone Screw Super Fastener ST	38mm
340-2040	4.0 Bone Screw Super Fastener ST	40mm
340-2042	4.0 Bone Screw Super Fastener ST	42mm
340-2044	4.0 Bone Screw Super Fastener ST	44mm
340-2046	4.0 Bone Screw Super Fastener ST	46mm
340-2048	4.0 Bone Screw Super Fastener ST	48mm
340-2050	4.0 Bone Screw Super Fastener ST	50mm
340-2052	4.0 Bone Screw Super Fastener ST	52mm
340-2054	4.0 Bone Screw Super Fastener ST	54mm
340-2056	4.0 Bone Screw Super Fastener ST	56mm
340-2058	4.0 Bone Screw Super Fastener ST	58mm
340-2060	4.0 Bone Screw Super Fastener ST	60mm
340-2062	4.0 Bone Screw Super Fastener ST	62mm
340-2064	4.0 Bone Screw Super Fastener ST	64mm
340-2066	4.0 Bone Screw Super Fastener ST	66mm
340-2068	4.0 Bone Screw Super Fastener ST	68mm
340-2070	4.0 Bone Screw Super Fastener ST	70mm
340-2072	4.0 Bone Screw Super Fastener ST	72mm
340-2074	4.0 Bone Screw Super Fastener ST	74mm
340-2076	4.0 Bone Screw Super Fastener ST	76mm
340-2078	4.0 Bone Screw Super Fastener ST	78mm
340-2080	4.0 Bone Screw Super Fastener ST	80mm
340-2085	4.0 Bone Screw Super Fastener ST	85mm
340-2090	4.0 Bone Screw Super Fastener ST	90mm

4.0 Super Fastener Instruments

Part Number	Description	Quantity
321-0001	4.0mm 2.5mm Drill Sleeve	1
321-0003	4.0mm Drill 110mm	2
321-0004	4.0mm Drill 200mm Calibrated	2
321-0005	Depth Gauge	1
321-0006	4.0 SF Bone Screw Tap, 145mm	1
110500	3.5mm Countersink	1
110702	2.5mm Hex Driver with Silicone Handle	1
110103	2.5mm Hex Driver w/ AO Adapter	1
110104	Long Hex Driver	1
110310	QC T-Handle	1
110013	2.5mm Drill Bit 110mm	2
110040	2.5mm Drill Bit 200mm	2
110205	Depth Gauge, 60mm	1
321-90001	4.0 Super Fastener Tray Outer	1
321-90002	4.0 Super Fastener Tray Caddy	1
321-90003	4.0 Super Fastener Tray Lid	1

Note

This document is intended exclusively for physicians and is not intended for laypersons. Information on the products and procedures contained in this document is of a general nature and does not represent and does not constitute medical advice or recommendations. Because this information does not purport to constitute any diagnostic or therapeutic statement with regard to any individual medical case, each patient must be examined and advised individually, and this document does not replace the need for such examination and/or advice.

All content herein is protected by copyright, trademarks and other intellectual property rights owned by or licensed to OsteoCentric Technologies, Inc. or one of its affiliates and must not be redistributed, duplicated, or disclosed, in whole or in part, without the express written consent of OsteoCentric Technologies, Inc.



75 West 300 N, Suite 150
Logan UT, 84321
Phone: 1-800-969-0639
info@osteocentric.com
osteocentric.com